

BAUSCH + LOMB

◀bi▶ Blade®



Experience the Dual-Port Difference

We are constantly seeking new ways to put innovation into your practice, and the Bi-Blade® vitrectomy cutter is an example of our dedication to providing a precise and efficient surgical experience.

STABILITY

15,000 CPM*

The Bi-Blade® vitrectomy cutter's dual-edge blade doubles the effective cut rate over single-edge vitrectomy cutters, allowing for faster procedures. With the consistent flow provided by the dual-port design and 100% open duty cycle, surgeons can feel confident in the ability to shave near a mobile retina, perform dissections and remove intraocular tissues.


*Effective cut rate.

EFFICIENCY

UP TO 230% FASTER FLOW RATE**

The dual-port greatly enhances efficiency due to its unique design that allows for cutting without port closure. As the port remains open, Bi-Blade® provides continuous aspiration at all times while cutting. We have reduced the need for a duty cycle while creating a consistent flow rate and continuous holding force.

**Compared to traditional single-port cutters in BSS on systems in 7500 CPM mode.

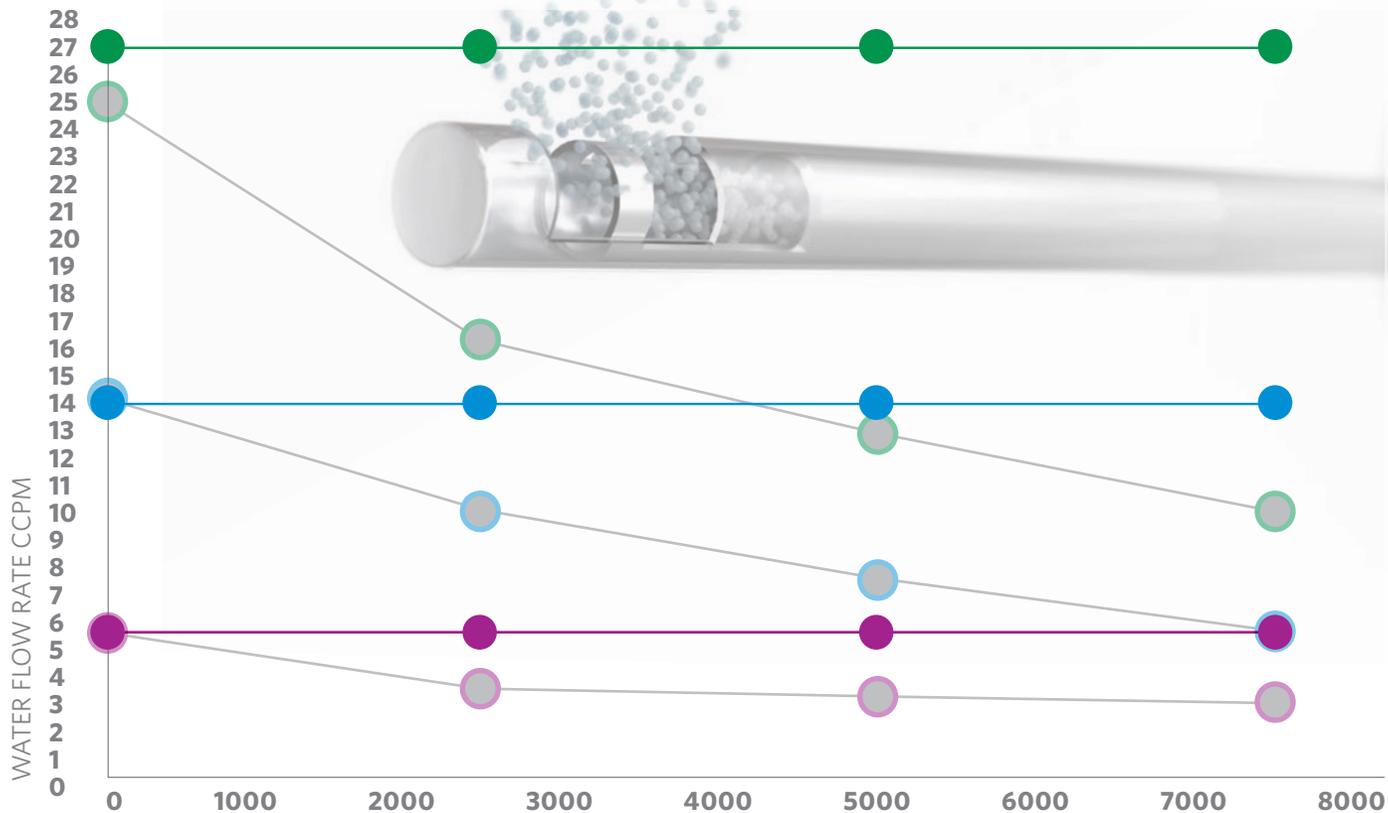

Stellaris Elite® CPM MODE

23ga Bi-Blade® at max flow in BSS — 23ga single-port cutter at max flow in BSS

25ga Bi-Blade® at max flow in BSS — 25ga single-port cutter at max flow in BSS

27ga Bi-Blade® at max flow in BSS — 27ga single-port cutter at max flow in BSS

◀bi▶ Blade®

Bi-Blade® Vitrectomy Cutter,
Stellaris Elite®

23G SE5623B
25G SE5625B
27G SE5627B



6/BOX UNLESS NOTED OTHERWISE

BAUSCH + LOMB

Bi-Blade® Vitrectomy Cutter Stellaris Elite® Packs

Posterior Mid-Field

23G SE5423MVB
25G SE5425MVB
27G SE5427MVB



Posterior Wide Field

23G SE5423WVB
25G SE5425WVB



Combined Mid-Field

23G SE5523MVB
25G SE5525MVB
27G SE5527MVB

Combined Wide Field

23G SE5523WVB
25G SE5525WVB



IMPORTANT SAFETY INFORMATION

Indications

The Bi-Blade® Vitrectomy Cutter is a pneumatically actuated device that supports aspiration and guillotine-style cutting functions to remove vitreous and/or other tissues from the eye during surgery. The Bi-Blade® Vitrectomy Cutter and the Bi-Blade® Vitrectomy Cutter Packs are intended for use in conjunction with the Stellaris Elite® or Stellaris PC® Vision Enhancement Systems, as specified.

Warnings

- Failure to follow priming instructions for the vitrectomy cutter and aspiration tubing may result in poor aspiration response and reduction in vitrectomy cutter's effectiveness.
- Do not re-sterilize and do not re-use.
- Only use the Bi-Blade® Vitrectomy Cutter and Bi-Blade® Vitrectomy Cutter Packs with the Stellaris Elite® or Stellaris PC® Vision Enhancement Systems, as specified.

This is not all you need to know. See the Bi-Blade® instructions for use for detailed directions, proper use, and full risk and safety information.

CAUTION: Federal (U.S.) Law restricts these devices to sale, by or on the order of a physician.

▶ Schedule a demo at [bauschretina.com](https://www.bauschretina.com)

BAUSCH + LOMB

Stellaris Elite is a registered trademark of Bausch & Lomb Incorporated or its affiliates. Bi-Blade is a registered trademark of Medical Instrument Development Laboratories, Inc. used under license. ©2021 Bausch & Lomb Incorporated or its affiliates. SRG.0028.USA.21