

enVista®

Hydrophobic Acrylic IOL

MONOFOCAL

All monofocals are
not created equal.
Make yours enVista®.
The clear choice for excellent
predictable outcomes.

ABERRATION-FREE
OPTIC

GLISTENING-FREE¹

CRAFTED FOR EXCEPTIONAL TOLERANCE TO
DECENTRATION & TILT^{2,3}

The TruSight® advanced optic (AO) with SureEdge® design is scratch-resistant and 25X the hardness of traditional hydrophobic acrylic IOLs.^{1*}

*Based on bench testing with a model cornea.

Aberration-free optics complement the natural depth of focus of the cornea and promotes excellent visual outcomes for your patients.²


StableFlex® technology supports efficient lens unfolding, ensuring excellent optic recovery.¹

enVista®
HYDROPHOBIC ACRYLIC IOL

BAUSCH + LOMB



MX60PL order number MXUEPLXXXX
MX60E / EE order number MX60UEXXXX / EEUXXXX

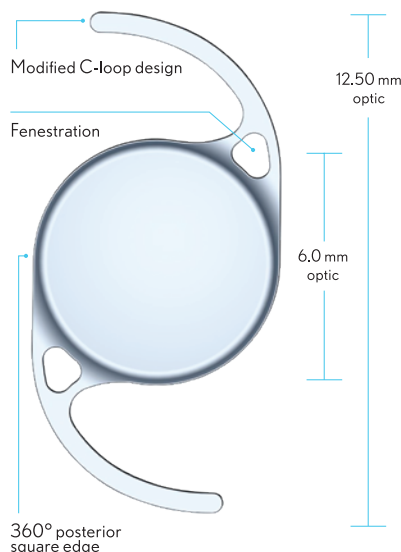

MODEL NUMBER	MX60PL (preload)	MX60E / EE (non-preload)
OPTIC DESIGN	One-piece Hydrophobic acrylic Aspheric, aberration-free, biconvex	One-piece Hydrophobic acrylic Aspheric, aberration-free, biconvex
OPTIC SIZE	6mm	6mm
LENGTH	12.5mm	12.5mm
HAPTICS	Modified C, fenestrated	Modified C, fenestrated
OPTICAL BIOMETRY	Optical A-constant* ACD Surgeon Factor	Optical A-constant* ACD Surgeon Factor
APPLANATION BIOMETRY	Applanation A-constant* ACD Surgeon Factor	Applanation A-constant* ACD Surgeon Factor
OTHER FEATURES	Glistening free Refractive index: 1.53 UV absorbing Sharp 360° square posterior edge	Glistening free Refractive index: 1.53 UV absorbing Sharp 360° square posterior edge
DIOPTER RANGE	0 to +10 D in 1.0-D increments +10 to +30 D in 0.5-D increments +30 to +34 D in 1.0-D increments	0 to +10 D in 1.0-D increments +10 to +30 D in 0.5-D increments +30 to +34 D in 1.0-D increments

* A-constant values are suggested as a guideline. Physicians should calculate lens power based on optimization of their experience and preference with IOL technology.

STORZ®
Ophthalmic Instruments

Find B+L IOL surgical equipment
online at www.StorzEye.com

Storz® BLIS Inserter System



FOR INSERTING LENS MODEL MX60E; +0 D to +34 D with X1 cartridge; +0 D to +34 D with X0 cartridge
RECOMMENDED INCISION SIZE BLIS-X1 2.4mm or less; BLIS-X0 2.2mm-2.6mm
TYPE OF ACTION Screw-type
COMMENTS Controlled delivery. Reusable. Sterilization required.

INJ100 Delivery System



FOR INSERTING LENS MODEL MX60PL
RECOMMENDED INCISION SIZE 2.2mm-2.6mm
TYPE OF ACTION Silicone tip push-type
COMMENTS Single-handed delivery. Disposable.

1. enVista directions for use.
2. Packer M. enVista hydrophobic acrylic intraocular lens: glistening free. Expert Review of Ophthalmology. 2015; 10:5,415-420.
3. Altmann GE, Nichamin LD, Lane SS, Pepose JS. Optical performance of 3 intraocular lens designs in the presence of decentration. J Cataract Refract Surg. 2005 Mar;31(3):574-85

Indications and Important Safety Information for enVista® IOL

INDICATIONS: Indicated for primary implantation for the visual correction of aphakia in adult patients in whom the cataractous lens has been removed. The lens is intended for placement in the capsular bag.

WARNINGS: Careful preoperative evaluation and sound clinical judgment should be used by the surgeon to decide the risk / benefit ratio before implanting a lens in a patient.

PRECAUTIONS: Do not resterilize this intraocular lens by any method. Do not use if the packaging is damaged or if there are signs of leakage. Do not store lenses at temperatures over 43°C (109°F) or lower than 0°C (32°F). Do not reuse the lens. Safety and effectiveness of the enVista® IOL have not been substantiated in patients with conditions and intraoperative complications as outlined in the enVista® IOL Directions for Use.

ADVERSE EVENTS: As with any surgical procedure, there is risk involved. Potential complications accompanying cataract or implant surgery may include, but are not limited to the following: corneal endothelial damage, infection (endophthalmitis), retinal detachment, vitritis, cystoid macular edema, corneal edema, pupillary block, cyclitic membrane, iris prolapse, hypopyon transient or persistent glaucoma, and secondary surgical intervention.

CAUTION: Federal law restricts this device to sale by or on the order of a physician.

ATTENTION: Reference the Directions for Use labeling for a complete listing of indications and important safety information.